


Air-Pro™ Series

Alternating Pressure Pad & Pump Overlay Mattress Systems

AIR-PRO™ SERIES - OVERLAY MATTRESS SYSTEM

The Air-Pro™ Plus & Elite Alternating Pump and Overlay Mattress Systems include a foldaway bed hanger bracket. This feature keeps the pump off the floor and reduces any accidental damage. Our easy access filter located on the bottom of the pump reduces any chance of dust particles from entering the pump compressor motor. We include a heavy duty, deluxe bubble style pad with head and foot tuck under flaps. The adjustability of the Elite allows this system to be used for heavier patients. Air-Pro™ meets California Technical Bulletin #117 for fire retardancy.

Shown with 4810 pad

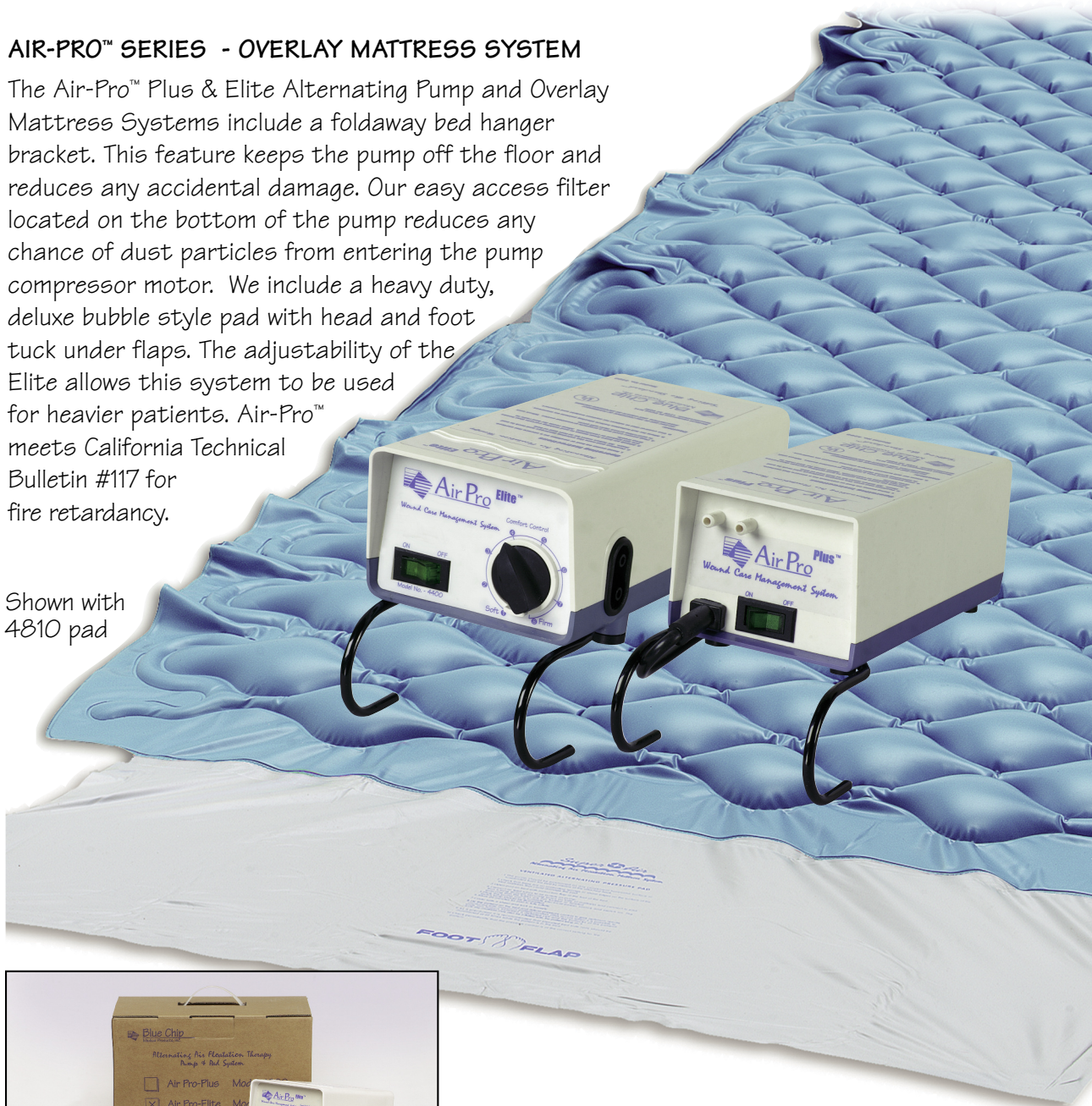

Model 4400 system shown

- Air-Pro™Plus Model 4200 with flaps
- Air-Pro™Plus Model 4200S without flaps
- Air-Pro™Elite Model 4400 with flaps
- Air-Pro™Elite Model 4400S without flaps


AIR-PRO™

SPECIFICATIONS

AIR-PRO™ SERIES - OVERLAY MATTRESS PAD

Material: Medical grade PVC with Anti-microbial, fire-retardant treatments.

Design: Bubble Style with more touch points for effective pressure redistribution.

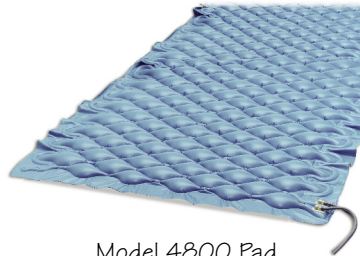


Model 4800 Pad without flaps



Model 4810 Pad with flaps

	Air-Pro™ Pad Solo	Air-Pro™ Pad Deluxe	Air-Pro™ Pad Laser	Air-Pro™ Pad Multi
Model:	4800	4810	4820	4830
Size:	35" x 79" x 2.5"	35" x 79" x 2.5"	35" x 79" x 2.5"	35" x 79" x 3.75"
End Flaps:	No	Yes	Yes	Yes
Ventilated:	No	No	36 Lasered Holes	36 Lasered Holes
Pump:	4201/4401	4201/4401	4401SS	4401SS
Heavy Duty:	Yes	Yes	Yes	Yes

MODEL 4200 System - 4201 Pump w/ 4810 Overlay Pad
MODEL 4200S System - 4201 Pump w/ 4800 Overlay Pad
MODEL 4201 Pump Only


Dimensions: 9L" x 4W" x 3.5H"
Weight: 4 lbs.
Electrical: 110V 60 hz
Listings: 
Cycle Time: 10 Minutes
Air Output: 5.6 Liters

*All models available in 220-240V; 50-60Hz



Model 4201

MODEL 4400 System - 4401 Pump w/ 4810 Overlay Pad
MODEL 4400S System - 4401 Pump w/ 4800 Overlay Pad
MODEL 4401 Pump Only
MODEL 4401SS Pump Only (to be used with Laser & Multi pad only)

Dimensions: 10L" x 5W" x 4H"
Weight: 5 lbs.
Electrical: 110V 60 hz
Listings: 
Cycle Time: 10 Minutes
Air Output: 6.5 Liters
Air Output (4401S): 8 Liters

*All models available in 220-240V; 50-60Hz



Model 4401



Model 4401SS
To be used with Laser & Multi pad only

PDAC Approved
Medicare E0181

WARRANTY

Pump: 1 year limited
Mattress: 6 months limited



Blue Chip

Medical Products, Inc.

7-11 Suffern Place, Suffern NY 10901

Tel: (845) 369-7535 • Fax: (845) 369-7633 • Toll Free: (800) 795-6115

www.bluechipmedical.com

"Setting the Standard"