

SEATING

GERI-GEL OVERLAY - MODEL 6200NS

Patient comfort and pressure redistribution is achieved for individuals who sit for extended periods of time during the day. Fits most standard recliners. Special sizes available for use in Dialysis and Chemotherapy recliners. Ideal for home style lift chairs. The three section overlay is constructed with an aqueous, highly viscous gel that is encapsulated in a high density foam shell. The overlay is covered with Vyvex®-III Stretch Cover and a non-skid Bottom. An adjustable safety strap secures the overlay for added safety. Cover material is anti-microbial and fire-retardant. Supports patients weights of 250 pounds. Bariatric and custom sizes available.

ZIPPERED COVER	CORE
Top: Vyvex® III - Multi-stretch knit urethane Low shear, fluid proof and stain proof	Top: Med-Flex® PPF® Center: Highly Viscous Gel Base: High Density Med-Flex® PPF®
Bottom: Non-skid tri-laminate vinyl	

Limited 12 Month Warranty

18" W x 70" L x 1.75" H



19" W x 70" L x 4" H
29" W x 70" L x 4" H

Custom sizes and weights available upon request.

GERI-GEL BARIATRIC OVERLAY - MODEL BCGG-6000BAR

ZIPPERED COVER	CORE
Top: Vyvex® III - Multi-stretch knit urethane Low shear, fluid proof and stain proof	Top: Gel-Flex® incorporated into Reflect® PPF® center section Base: High Density Med-Flex® PPF®
Bottom: Non-skid tri-laminate vinyl	

Standard size supports patient's weight up to 450 lbs.
Can be manufactured to support patients weights up to 800 lbs.

Limited 12 Month Warranty

CHAIR-AIR® MODEL 9700GR GERIATRIC RECLINER MATTRESS & PUMP SYSTEM

This unique alternating pressure mattress system allows for patients who are at high risk for pressure sores to be placed in a variety of mobile devices. Patients are no longer confined to bed for the whole day due to risks of sitting without proper pressure redistribution. Chair-Air® systems can be integrated into a geriatric chair, recliner or standard high-back wheelchair. Twenty bladders gently alternate, providing low interface pressures in addition to comfortably supporting the patient. The Chair-Air® 9700GR system includes the 9701 alternating pump. Chair-Air® meets California Technical Bulletin #117 for fire retardancy.

ZIPPERED COVER	CORE 3.25"
Top: Vyvex® III - Multi-stretch knit urethane Low shear, fluid proof and stain proof	Bladders: 20 - 3.25" air cells

Limited 12 Month Warranty



Model 9701




25" x 45" x 3.25"

CHAIR-AIR® MODEL 9700C & 9700CAF ALTERNATING AIR WHEELCHAIR CUSHION AND PUMP



9700C

Shown in cover

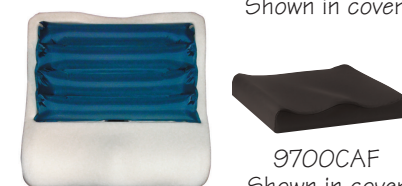

9700CAF

Shown in cover

The Chair-Air® 9700C has six individual bladders that alternate to provide the patient a superior pressure redistribution seating surface. The zippered non-skid stretch-knit cover is anti-microbial and fire-retardant. The Chair-Air® 9700C includes the Chair-Air® 9701 alternating pump.

ZIPPERED COVER	CORE 9700C 3.25", 9700CAF 3"
Top: Low shear stretch knit anti-microbial nylon	Bladders 9700C: 6 - 3.25" air cells Bladders 9700CAF: 4 - 2.5" air cells

Limited 12 Month Warranty

The Chair-Air® 9700CAF has four individual bladders that alternate to provide the patient a superior pressure redistribution seating surface while the unique molded foam design, with slight abductorommel, stabilizes the hips and keeps the patient in proper mid-line alignment. The zippered non-skid stretch-knit cover is anti-microbial and fire-retardant. The Chair-Air® 9700CAF includes the Chair-Air® 9701 alternating pump.

WHEELCHAIR POSITIONING BELTS



Safety is a must for all patients sitting in a wheelchair. We offer three styles of belts that are easy to attach to any standard wheelchair frame.

Auto Push Button: This buckle style offers a spring loaded "pop-out" connector. The center release mechanism separates with minimum effort.

Side Release: This injection molded Delrin "side release" buckle provides excellent durability with light weight structure. Just squeeze the sides for quick release.

Hook and Closure: The remarkable sheer strength of hook and loop closure plus its "Cinch" features makes this seat belt one of our best. It is easy to use, lightweight and boasts a simple on and off quality.

MODEL	DESCRIPTION
550-AB	Auto Push Button
570-SR	Side Release
580-VC	Hook and Closure
CUSTOM SIZES AVAILABLE	

ANTI-THRUST

ANTI-THRUST - MODEL 402

Provides the proper positioning for patients who hyper-extend by thrusting forward. Cushion base and top portion is a high density foam and is offered with a pressure redistribution foam top. Gel can be added for additional pressure redistribution. Covered with medical grade vinyl, non-skid bottom with safety strap.

ZIPPERED COVER	CORE
Cover: Medical Grade Vinyl Option: Vyvex® III - Multi-stretch knit urethane Low shear, fluid & stain proof, incontinent	Core: Med-Flex® PPF® high density foam Option: Viscous Gel Model 4020

Limited 12 Month Warranty

WEDGE- MODEL 608

Wedge cushions significantly reduce the patient's ability to hyperextend out of the wheelchair. The cushion is a high density foam that increases the seat to back angle, just enough to maintain proper sitting posture. This cushion is available with gel for pressure redistribution and immersion. The standard cover is a gray vinyl. Also available in black Vyvex®-III stretch top with non-skid bottom. All cover materials are anti-microbial and meet California Technical Bulletin #117 for fire retardancy. An optional solid seat insert is available as well as safety strap and buckle. Supports 250 lbs.

ZIPPERED COVER	CORE 2"- 1", 3"- 1", 4"- 2"
Cover: Medical Grade Vinyl Option: Vyvex® III - Multi-stretch knit urethane Low shear, fluid & stain proof, incontinent	Top: High Density Med-Flex® PPF® Option: Viscous Gel Model 6080
Option: Non-skid tri-laminate vinyl Safety strap and buckle	

Limited 12 Month Warranty

LEVEL-R SERIES MODEL 705

Internal rotation of the hips are reduced as well as a positive midline sitting position is achieved with this cushion. The "Level-R" reduces the hammock effect usually found with a wheelchair sling seat. Available in a wedge design to reduce patients from sliding out of the wheelchair. The "Level-R" is covered in a fire retardant, anti-microbial black vinyl cover. All cover materials meet California Technical Bulletin #117 for fire retardancy. Supports 250 lbs.



ZIPPERED COVER	CORE 2" TAPER TO 1"
Cover: Medical Grade Vinyl Option: Vyvex® III - Multi-stretch knit urethane Low shear, fluid & stain proof, incontinent	Core: Med-Flex® PPF® super high density foam
Option: Safety strap and buckle	


Limited 12 Month Warranty

LEVEL-R WEDGE MODEL 805



ZIPPERED COVER	CORE 3" TAPER TO 1.5", 4" TAPER TO 1.5"
Cover: Medical Grade Vinyl Option: Vyvex® III - Multi-stretch knit urethane Low shear, fluid & stain proof, incontinent	Core: Med-Flex® PPF® super high density foam
Option: Safety strap and buckle	

Limited 12 Month Warranty

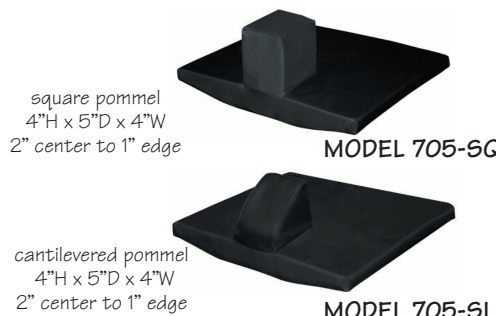


squareommel
4"H x 5"D x 4"W
2" center to 1" edge

MODEL 705-SQ

cantileveredommel
4"H x 5"D x 4"W
2" center to 1" edge

MODEL 705-SL

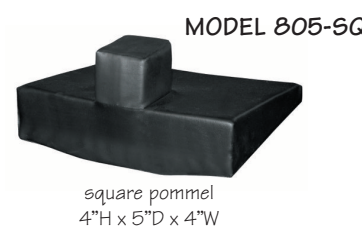

MODEL 805-SQ

squareommel
4"H x 5"D x 4"W



MODEL 805-SL

cantileveredommel
4"H x 5"D x 4"W



Wheelchair Backs



Wheelchair Cushions



Positioning Cushions



Anti-Thrust Cushions



Recliner Cushions

Seating & Positioning


SIZE	16" x 16"	17" x 17"	18" x 16"	20" x 16"	22" x 16"	24" x 16"	18" x 18"	20" x 18"	22" x 18"	24" x 18"	20" x 20"	24" x 20"
SUFFIX	NA		A	FA	VFA	VVFA	C	CC	CCFA	CCVFA	D	DD

SEATING & POSITIONING

GEL-PRO® POSITION PLUS® - MODEL 6250

Skin Protection and Position Cushion

The Gel Pro® Position Plus® Cushion provides the most advanced pressure redistribution, therapeutic support and positioning available. Through proper immersion of the bony prominences the Gel Pro® Position Plus® cushion allows therapists to offer correct posture, proper spinal alignment, trunk stability and all the pressure redistribution characteristics required for a seating application. All cover materials are anti-microbial and meet California Technical Bulletin #117 for fire retardancy. An optional solid seat insert is available. Supports 250 lbs.

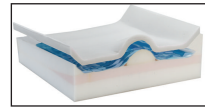

ZIPPERED COVER		CORE 4"	
Top: Vyvex® III - Multi-stretch knit urethane Low shear, fluid proof and stain proof	Top: Viscotec® memory foam	Center: Gel and PPF® high density foam	
Option: Vyvex® III - black	Base: Med-Flex® PPF® high density foam		
Bottom: Non-skid tri-laminate vinyl			
PDAC APPROVED Medicare E2607 - less than 22" • Medicare E2608 - 22"W OR GREATER			
Limited 18 Month Warranty			



PRO POSITIONER® - MODEL 6240

Skin Protection and Position Cushion




Pro Positioner® Cushion provides advanced pressure redistribution, therapeutic support and positioning. Through proper immersion of the bony prominences the Pro Positioner® cushion allows therapists to offer correct posture, proper spinal alignment, trunk stability and all the pressure redistribution characteristics required for a seating application. All cover materials are anti-microbial and meet California Technical Bulletin #117 for fire retardancy. An optional solid seat insert is available. Supports 250 lbs.

ZIPPERED COVER		CORE 4"	
Top: Vyvex® III - Multi-stretch knit urethane Low shear, fluid proof and stain proof	Top: Viscotec® memory foam	Center: PPF® high density foam	
Option: Vyvex® III - black	Base: Med-Flex® PPF® high density foam		
Bottom: Non-skid tri-laminate vinyl			
PDAC APPROVED Medicare E2607 - less than 22" • Medicare E2608 - 22"W OR GREATER			
Limited 18 Month Warranty			

GEL-PRO® ELITE - MODEL 6204

Skin Protection

Our segmented gel bladder contains an aqueous, highly viscous gel which provides for pressure redistribution, proper immersion and envelopment around the bony prominences. The bladder is surrounded by a high density foam shell that is covered using Vyvex-1®, a fluid-proof, low shear nylon top with a medical grade vinyl bottom. Supports 250 lbs.

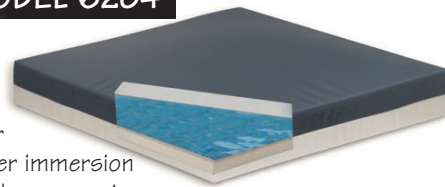

ZIPPERED COVER		CORE 2.75"	
Top: Vyvex® I - Low shear, fluid proof anti-microbial nylon	Top: Med-Flex® PPF®	Center: Ultra soft viscous gel bladder	
Bottom: Medical grade vinyl	Base: Med-Flex® PPF® high density foam		
PDAC APPROVED Medicare E2603 - less than 22" • Medicare E2604 - 22"W OR GREATER			
Limited 18 Month Warranty			

COMFORT CARE® - MODEL 602-3MC

General Use

This high density foam cushion is covered in Vyvex I®, a low shear fluid-proof, anti-microbial nylon cover. Supports 250 lbs.



ZIPPERED COVER		CORE 3"	
COVER: Vyvex® I - Low shear, fluid proof anti-microbial nylon	CORE: Med-Flex® PPF® high density foam		
Option: Vyvex® III, Multi-stretch			
PDAC APPROVED Medicare E2601 - less than 22" • Medicare E2602 - 22"W OR GREATER			
Limited 18 Month Warranty			

PRODUCT ORDERING GUIDE Add suffix to model number (ex. Model 6204 + Size 16" x 16" = 6204N)

SIZE	16" x 16"	17" x 17"	18" x 16"	20" x 16"	22" x 16"	24" x 16"	18" x 18"	20" x 18"	22" x 18"	24" x 18"	20" x 20"	24" x 20"
SUFFIX	NA		A	FA	VFA	VVFA	C	CC	CCFA	CCVFA	D	DD

SEATING & POSITIONING

AMARA® Component based system allows for customization.

Amara® is a component based system which offers a variety of positioning bases in conjunction with a selection of pressure redistribution inserts. Amara's innovative deep well design is incorporated into the cushion's high density Medflex® or Viscotec® visco elastic foam base under the ischials and coccyx. The rehabilitation specialists choice of pressure redistribution Viscotec® visco elastic foam, Soft-Flo® viscous gel or Gel-Flex® solid gel inserts redistribute pressure and reduce shearing forces.

ZIPPERED COVER: Top - Vyvex® III - Urethane coated, multi-stretch low shear, fluid & stain proof, incontinent Bottom - Non-Skid



AMARA® 100 - GENERAL USE

Patients with low to moderate risk of skin breakdown and who require minimal postural support. Bariatric sizes available.

BASE: Med-Flex® PPF®

MODEL AMR-100	PDAC APPROVED Medicare E2601 - less than 22" Medicare E2602 - 22"W OR GREATER
	Limited 18 Month Warranty

AMARA® 300 - SKIN PROTECTION

For patients with moderate to high risk of skin breakdown and require minimal postural support. Bariatric sizes available.

BASE: MedFlex® PPF®

INSERT: Viscotec®

MODEL AMR-300	PDAC APPROVED Medicare E2603 - less than 22" Medicare E2604 - 22"W OR GREATER
	Limited 18 Month Warranty

AMARA® 450 - SKIN PROTECTION/POSITIONING

Patients with moderate to high skin breakdown and who require significant postural support.

BASE: Med-Flex® PPF®

INSERT: Soft-Flo®

MODEL AMR-450	PDAC APPROVED Medicare E2607 - less than 22" Medicare E2608 - 22"W OR GREATER
	Limited 18 Month Warranty

AMARA THERAPEUTIC WHEELCHAIR BACK SYSTEM

Amara® Therapeutic Wheelchair Back System improves the comfort and support of your existing wheelchair back by eliminating the hammock effect. Amara provides proper patient positioning to maintain spinal and shoulder alignment and superior pressure redistribution to prevent skin breakdown. Amara® comes with either a Dual-Zone design (BK-1), Fixed Back (BK-2) or Adjustable Lumbar Support (BK-3) which is easy to adjust. Amara® combines the unmatched support of Med-Flex® foam with the therapeutic pressure redistribution qualities of Viscotec® foam.

PDAC APPROVED
Medicare E2611 - less than 22"
Medicare E2612 - 22"W OR GREATER



PRODUCT ORDERING GUIDE Add suffix to model number (ex. Model AMR-450 + Size 16" x 16" = AMR-450-1616)


SIZE	16" x 16"	18" x 16"	16" x 18"	18" x 18"	18" x 20"	20" x 16"	20" x 18"	20" x 20"	22" x 18"	22" x 20"	24" x 18"	24" x 20"	26" x 20"	28" x 20"	30" x 20"
SUFFIX	1616	1816	1618	1818	1820	2016	2018	2020	2218	2220	2418	2420	2620	2820	3020

*AMARA® back available sizes in grey

SEATING & POSITIONING

SELF ADJUSTING GRZ® - MODEL GRZ-WC

3 independent air cell chambers

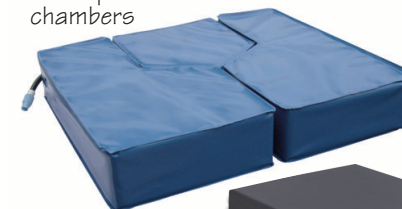

The GRZ® is dynamic adjustable air bladder system consisting of 3 individual air cell chambers which redistribute pressure and adjust to body weight and position. Ideal for patients with moderate to high risk of skin breakdown who require significant postural support. The GRZ® allows the therapist to offer correct posture, proper alignment and trunk stability as well as adjustments for pelvic obliquities. Does not require pump for adjustment. Adjustable strap secures cushion to wheelchair. Supports up to 300lbs.

PDAC APPROVED
Medicare K0734 - less than 22"
Medicare K0735 - 22"W OR GREATER

Limited 18 Month Warranty

ZIPPERED COVER		CORE 3"	
COVER: Vyvex® III - Multi-stretch knit urethane low shear, fluid & stain proof, incontinent		CORE: 3 individual air chambers	

COCCYX SEAT CUSHION - MODEL 603

The coccyx cutout redistributes pressure at the sacral area while offering patient comfort and even distribution of weight. Medium density fire-retardant foam is covered in either a black or navy poly/cotton or gray medical grade vinyl. An optional solid seat insert is available. Supports 250 lbs.



SEALED COVER		CORE 2", 3" & 4"	
COVER: Poly/COTTON	CORE: Medium Density PPF® FOAM		
Option: Medical grade gray vinyl			
Limited 12 Month Warranty			

VISCOTEC™ FOAM CUSHION - MODEL 604

Viscotec™ foam conforms to the patient's body by surrounding bony prominences, reducing interface pressures. This heat activated foam self-adjusts to the seating surface to maximize comfort and pressure redistribution. Viscotec™ foam is covered using either a 4-way stretch-knit or Vyvex III® low shear, fluid-proof four-way stretch material. Vyvex III® and foam are fire-retardant and anti-microbial. Navy stretch-knit cover is shipped unless otherwise specified.

FLAP-FOLD COVER		CORE 2", 3" & 4" (specify)	
COVER: STRETCH KNIT NAVY	CORE: VISCOTEC™ FOAM		
Option: Vyvex® III zippered Multi-stretch knit			
Limited 12 Month Warranty			




SOLID SEAT INSERT - MODEL SSI-600 SOLID BACK INSERT - MODEL SSI-BK-600

Solid seats and backs are used to properly align the body as well as to eliminate the sagging effect found in standard sling seat wheelchair upholstery. Safety straps secure the solid back to the wheelchair. The cushion uses medium density foam and is covered using waterproof black leatherette cover.

COVER		CORE 1.25"	
COVER: LEATHERETTE	CORE: Med-Flex® PPF® FOAM	BASE: Wooden Solid SEAT INSERT	
Limited 12 Month Warranty			

SOLID SEAT AMPUTEE CUSHION - MODEL SSA-600

This versatile cushion utilizes Med-Flex® PPF® foam over a wooden solid seat base. The low shear nylon cover is anti-microbial and fluid-proof to protect the inner core. Supports 250 lbs. Left or right 6" stump extension. Safety strap secures insert to the wheelchair.

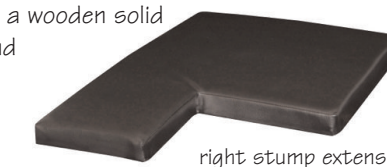

COVER		CORE 1.75"	
COVER: LEATHERETTE	TOP: Med-Flex® PPF® FOAM	BASE: Wooden Solid SEAT INSERT	
Option: Vyvex® III - Multi-stretch knit urethane Low shear, fluid & stain proof, incontinent Viscotec foam STUMP EXTENSION			
Limited 12 Month Warranty			



left stump extension

PRODUCT ORDERING GUIDE Add suffix to model number (ex. Model 6204 + Size 16" x 16" = 6204N)

SIZE	16" x 16"	17" x 17"	18" x 16"	20" x 16"	22" x 16"	24" x 16"	18" x 18"	20" x 18"	22" x 18"	24" x 18"	20" x 20"	24" x 20"
SUFFIX	NA		A	FA	VFA	VVFA	C	CC	CCFA	CCVFA	D	DD