

Gel-Pro® Series

Gel-Pro® Position^{PLUS}®

GEL-PRO® POSITION ^{PLUS}®

The Gel Pro® Position ^{PLUS}® Cushion provides the most advanced pressure relieving therapeutic support and positioning available. The Gel Pro® Position ^{PLUS}® Cushion allows therapists to offer correct posture, proper spinal alignment, trunk stability and all the pressure relieving characteristics required for a seating application. Available in low profile design.

The Vyvex III® top cover is a low shear, stretch urethane, which is fluid and stain proof. The bottom cover is a non-skid tri-laminate vinyl.



Model 6250

The front abduction wedge maintains proper midline alignment and reduces internal rotation of the hips in addition to preventing sliding and shearing forces.

The cushion's top layer is Viscotec®, a heat activated memory foam which surrounds and supports the bony prominences while reducing interface pressure.

The middle layer is VPF®, a high density pressure relieving foam, encapsulating an ultra soft viscous gel bladder further reducing interface pressure.

A Medflex® foam base provides the necessary foundation to ensure that the patient maintains the proper stability, support and correct spinal alignment.



All foams and cover materials are anti-microbial and meet California Technical Bulletin #117 for fire retardancy.

GEL-PRO®

SPECIFICATIONS

GEL-PRO® POSITION PLUS®

Zippered Cover:


Top: Vyvex III® - Four-Way Stretch Knit Urethane
Low shear, fluid proof and stain proof

Bottom: Non-skid tri-laminate vinyl



Core:

Top: Viscotec® Memory Foam
Center: VPF® High Density Foam
Ultra-Soft Viscous Gel Bladder
Base: Medflex® VPF® High Density Foam

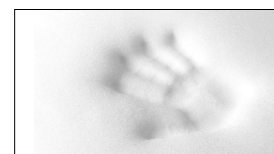

Heat Activated Viscotec®
Memory Foam

GEL-PRO® POSITION PLUS®

Model 6250NA 16" x 16"

Model 6250 17" x 17"

Model 6250A 18" x 16"

Model 6250FA 20" x 16"

Supports 250lbs.

Custom and Bariatric Sizes Available
Please state patients' weight when ordering

WARRANTY

Cushion: 1 year limited



7-11 Suffern Place, Suffern NY 10901

Tel: (845) 369-7535 • Fax: (845) 369-7633 • Toll Free: (800) 795-6115

www.bluechipmedical.com

"Setting the Standard"

Alternating Air-Floatation Mattresses

Alternating Pressure Pads & Pumps

Low Air-Loss Mattress Systems

Lateral Rotation Mattress Systems

Gel Mattresses

Gel Mattress Overlays

Static-Air Mattress Overlays

Water Mattress Overlays

Wheelchair Cushions

Patient Lift Slings

Side Rail Bumpers

Cart Covers

OR Replacement Table Mattresses

Replacement Stretcher Mattresses

Cubicle Curtains

Hospital Pillows