Relief-Care® Pro

Therapeutic Dual-Zone Foam Mattress with SCT" Technology

The Relief-Care® Pro Mattress System uses advanced Surface Cut Technology (SCT^M) to create a clinically effective support surface. SCT^M enables Blue Chip to produce a therapeutic dual-zone PPF® Pressurized Polymer Foam support surface with varying features and properties. This engineering advancement provides proper pressure redistribution that offers the patient a clinically effective mattress with low interface pressures. This design absorbs and redistributes patient's weight to provide greater comfort, maximum support and proper pressure redistribution around bony prominences.




The Relief-Care® Pro is covered with a Vyvex-I® top cover and a medical grade vinyl bottom cover. Vyvex-I® is a low shear, fluid proof nylon material which is anti-microbial, fire-retardant, stain resistant and washable. The foam core and cover meet Federal Flammability Standard 16 CFR 1633 for fire retardancy. Mattress can be manufactured to meet California Technical Bulletin #129 and Boston Fire Code 1X-11 when specified.

"SCT™ technology creates an advanced pressure redistribution support surface that provides unlimited comfort zones.

| MODEL | SIZE | THICKNESS |
|-----------|-----------|-----------|
| BCHTM 860 | 36" x 76" | 5.5" |
| BCHTM 880 | 36" x 80" | 5.5" |
| BCHTM 884 | 36" x 84" | 5.5" |
| BCHTM 842 | 42" x 80" | 7" |

CUSTOM AND BARIATRIC SIZES AVAILABLE WARRANTY: 12 MONTHS

PDAC Approved Medicare E0184

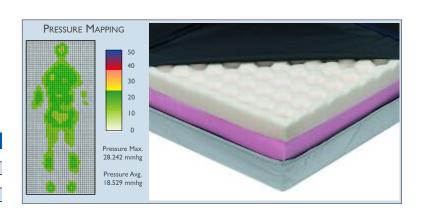




Home-Care Elite®

Homecare & Institutional High-Density Foam Mattress

Designed for use in Nursing Homes, Hospitals, and Institutions. The inner core is a high density, fire-retardant polyurethane foam that is covered with an anti-microbial, fire-retardant material meeting Federal Flammability Standard 16 CFR 1633. The mattress is closed with a self-ventilating zipper, and can be special ordered for use in prisons and psychiatric areas utilizing an inverted lock-stitch closure.

Single layer of high-density foam encased in a zippered medical grade vinyl cover providing superior patient comfort.

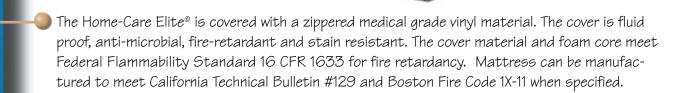



Durability, resiliency and patient comfort make the Home-Care Elite® the perfect mattress for homecare and institutional applications.

| MODEL | SIZE | THICKNESS |
|----------|-----------|-----------|
| BCHM 760 | 36" x 76" | 5.5" |
| BCHM 780 | 36" x 80" | 5.5" |
| BCHM 784 | 36" x 84" | 5.5" |
| BCHM 742 | 42" x 80" | 7" |

CUSTOM AND BARIATRIC SIZES AVAILABLE

WARRANTY: 12 MONTHS

