


GRZ™

Non-Powered Pressure Redistribution Mattress Replacement System

GRZ™

Blue Chip Medical's innovative GRZ™ helps to provide proper pressure redistribution to reduce interface pressures by dynamically adjusting to body weight and movement. The GRZ™ is a non-powered pressure redistribution mattress with nine individual dynamic self-adjusting air bladders inside a high resiliency foam shell with heel section. GRZ's nine zones redistribute pressure where it is needed most and supports patient weights up to 1000 lbs. with bariatric models.

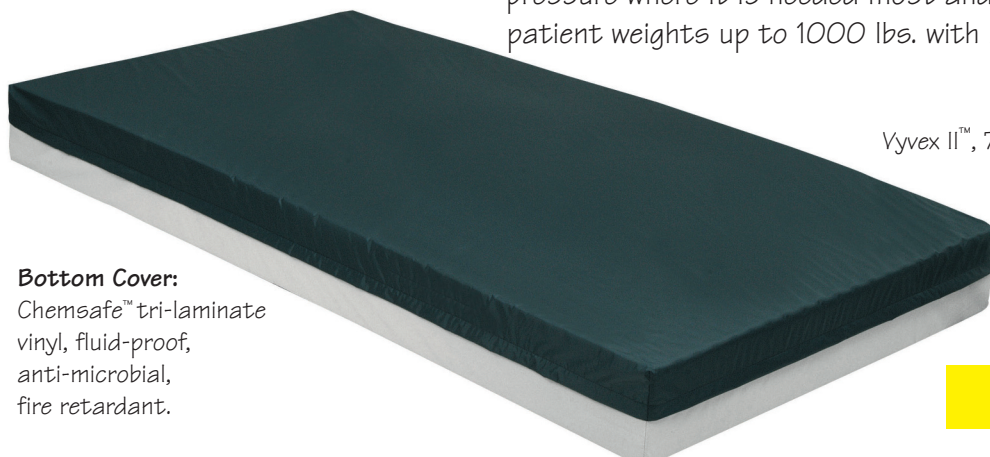


Top Cover:

Vyvox II™, 70 denier nylon, butyl backed, fluid-proof, anti-microbial, fire retardant.

Bottom Cover:

Chemsafe™ tri-laminate vinyl, fluid-proof, anti-microbial, fire retardant.


PDAC Approved
Medicare EO277

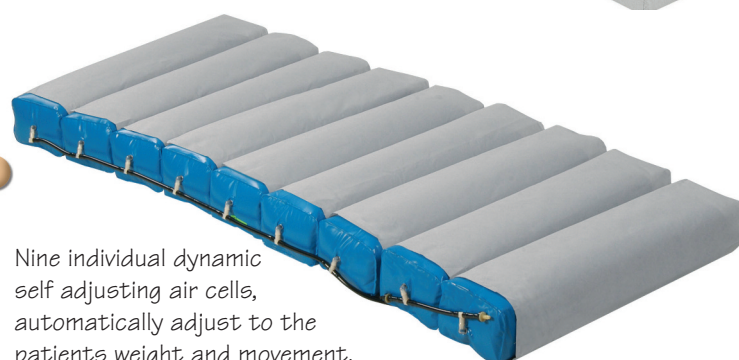

Nine individual dynamic self adjusting air cells, automatically adjust to the patients weight and movement. Urethane/vinyl blend/reforming foam agent, air circulating bladders. Enclosed sleeve for secure positioning.

Barrier Lok™ Infection Control System

Blue Chip's exclusive Barrier Lok™ feature virtually seals and protects the inner foam core of the mattress by preventing the leaching of fluids from penetrating the perimeter seam of the mattress. The Barrier Lok™ ensures that infection control standards are maintained, thereby providing investment protection on all GRZ™ mattress systems.

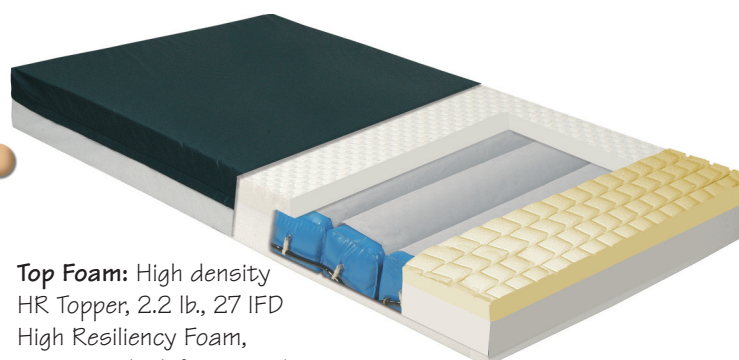

Top Foam: High density HR Topper, 2.2 lb., 27 IFD High Resiliency Foam, anti-microbial, fire retardant. Heel Wrap™ VFPT™ foam for heel protection.



Optional Pump

Alternating air pressure pump

 **Blue Chip**
Medical Products, Inc.

"Setting the Standard"

SPECIFICATIONS

MATTRESS ONLY

Dimensions:	36" x 76" x 7"
Weight:	48 lbs. (packaged)
Air Cells:	9 Self-adjusting air cells
Shell:	High resilience foam with Honeycomb SCT™ Heel Wrap™ VPF™ Foam

COVER

Top Cover:	Vyvetex II™, 70 denier nylon, butyl backed, fluid-proof, anti-microbial, fire retardant.
Bottom Cover:	Chemsafe™ tri-laminate vinyl, fluid-proof, anti-microbial, fire retardant.

Material meets California Technical Bulletin #117.

GRZ™ ORDERING INFORMATION

Model: GRZ 3676	Mattress 36" x 76" x 7"	Supports weight to 550
Model: GRZ 3680	Mattress 36" x 80" x 7"	Supports weight to 550
Model: GRZ 3684	Mattress 36" x 84" x 7"	Supports weight to 550
Model: GRZ 4280	Mattress 42" x 80" x 7"	Supports weight to 650
Model: GRZ 4880	Mattress 48" x 80" x 7"	Supports weight to 1000
Model: GRZ 5480	Mattress 54" x 80" x 7"	Supports weight to 1000
Model: GRZ 6080	Mattress 60" x 80" x 7"	Supports weight to 1000
Model: GRZ 4610	Optional Alternating Air Pressure Pump	

FLAMMABILITY INFORMATION

All GRZ™ Mattress Systems meet Federal Flammability Standard 16 CFR 1633. All mattresses can be ordered with a special upgraded fire barrier protection to meet flame retardancy specifications. This includes California Technical Bulletin #121, #129, 603 and City of Boston 1X-11.

WARRANTY

Mattress:	5 years limited
Cover:	6 months limited





7-11 Suffern Place, Suffern NY 10901

Tel: (845) 369-7535 • Fax: (845) 369-7633 • Toll Free: (800) 795-6115
www.bluechipmedical.com

"Setting the Standard"

©Blue Chip Medical Products, Inc. 2014

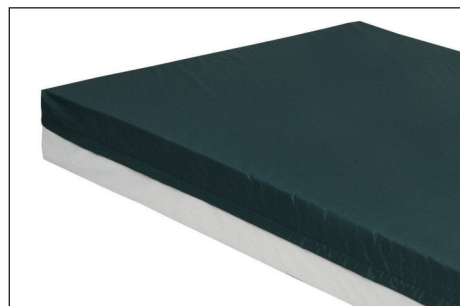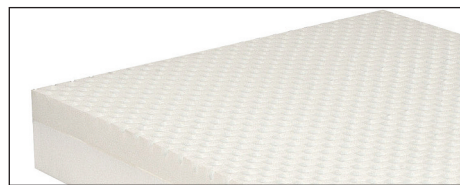

Alternating Air-Floatation Mattresses
Alternating Pressure Pads & Pumps
Low Air-Loss Mattress Systems
Lateral Rotation Mattress Systems
Gel Mattresses
Gel Mattress Overlays
Static-Air Mattress Overlays
Water Mattress Overlays
Wheelchair Cushions
Patient Lift Slings
Side Rail Bumpers
Cart Covers
OR Replacement Table Mattresses
Replacement Stretcher Mattresses
Cubicle Curtains
Hospital Pillows