


Support Surfaces

Pressure Redistribution Mattresses




Alternating Pressure Systems
with Low Air Loss

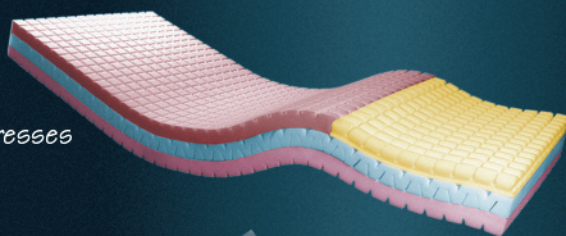

Lateral Rotation Systems
with Low Air Loss




Convertible Mattress Systems

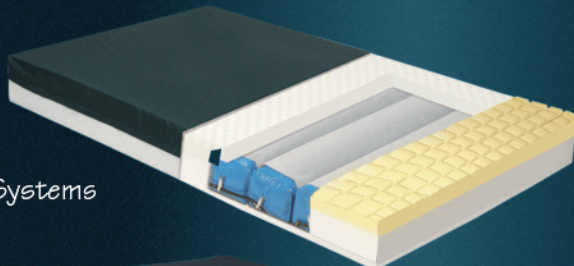

Homecare & Institutional Mattresses



Alternating Pressure Pads & Pumps

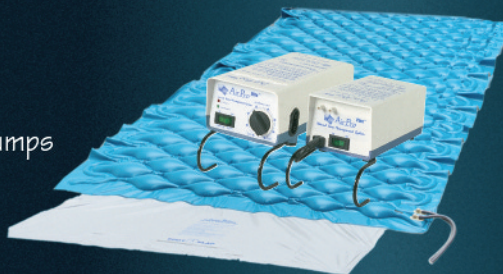

7-11 Suffern Place, Suffern NY 10901
Toll Free: (800) 795-6115 • Fax: (845) 369-7633
www.bluechipmedical.com

"Setting the Standard"

ALTERNATING PRESSURE WITH LOW AIR LOSS



Model 4600

APOLLO 3-PORT™

The Apollo 3-Port™ Alternating Pressure Low Air Loss Mattress System combines the most advanced technology and clinically effective features in a portable system. Provides pressure redistribution for the prevention and treatment of Stages I-IV pressure ulcers. Model 4600SRS with Safety Rail System.

Model:	4600
Supports:	350 lbs
Size:	35"W x 80"L x 8"H
Bladders:	18 - 8" air cells, ventilated, low air loss
Fire Code:	Cal. Tech Bulletin #117
Cover:	***Vyvex III™
Warranty:	1 Year Limited
PDAC:	E0277



Model 9600

SUPREME AIR™

The Supreme Air™ Alternating Pressure Low Air Loss Mattress System combines the most advanced technology with clinically effective static and alternating features in a portable system. Provides pressure redistribution for the prevention and treatment of Stages I-IV pressure ulcers. Includes low pressure alert and audible power outage alarm. Model 9600SRS with Safety Rail System.

Model:	9600
Supports:	350 lbs
Size:	35"W x 80"L x 8"H
Bladders:	18 - 8" air cells, ventilated, low air loss
Fire Code:	Cal. Tech Bulletin #117
Cover:	***Vyvex III™
Warranty:	1 Year Limited
PDAC:	E0277



Ventilated
Foot Section
9600 & 4600



Optional Digital Pump

POWER PRO™ ELITE

Combines two clinically effective therapies: Alternating Pressure with Low Air Loss in a portable blower based system. Provides proper pressure redistribution for the prevention and treatment of Stages I-IV pressure ulcers. The new Power-Pro VW Elite bariatric model fits 35", 39", 42" and 48" bed frames. An optional pulsating low air loss therapy feature is available. Model 9500SRS with Safety Rail System.



Model 9500

Model:	9500
Supports:	650 lbs. - 1,000 lbs.
Size:	36"W x 80"L x 9"H
	Available in 39"W, 42"W, 48"W
Bladders:	18 - 9" air cells, ventilated, low air loss
Fire Code:	Cal. Tech Bulletin #117
Cover:	***Vyvex III™
Warranty:	1 Year Limited
PDAC:	E0277

TRADEWIND™ SERIES

The Tradewind™ series offers therapists a choice of three dual mode pumps. All Tradewind™ pumps feature Alternating Pressure and Static modes as standard. The Alternating Pressure Mode provides the therapeutic benefit of reactive hyperemia and redistributes pressure well below 32mmHg, while the Static Mode provides gentle low air loss with constant pressure redistribution and support. This allows nursing to perform patient care.



TRW-1000
Mechanical pump

TRW-3000
Digital Pump


TRW-ATS
Digital Pump 750 lbs

Model:	TRW-1000 TRW-3000 TRW-ATS
Supports:	350 lbs. - 750 lbs.
Size:	35"W x 80"L x 8"H avail. in 42"W, 48"W
Bladders:	17 - 8" air cells, ventilated, low air loss
Fire Code:	Cal. Tech Bulletin #117
Cover:	***Vyvex III™
Warranty:	1 Year Limited
PDAC:	E0277

AIR SYSTEMS

GRZ™ - NON-POWERED CONVERTIBLE SYSTEM

The GRZ™ is a non-powered, pressure redistribution, self-adjusting, convertible mattress. GRZ™ helps to provide proper pressure redistribution to reduce interface pressures by dynamically adjusting to body weight and movement. The GRZ™ air cells are encased in a Medflex™ pressure redistribution PPF™ foam shell to redistribute pressure where it is needed most. Optional pump converts GRZ™ to alternating air therapy. Maintenance Free.

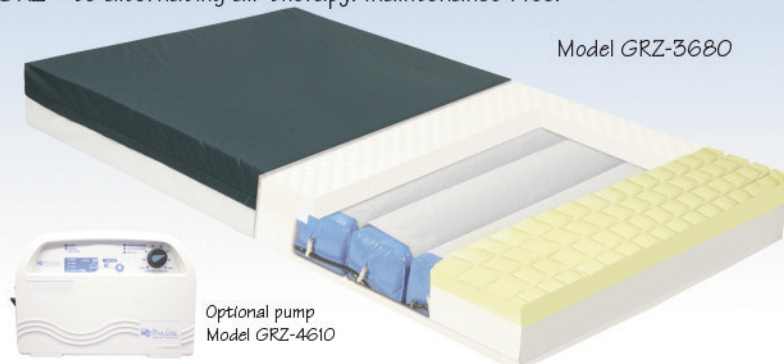

Model GRZ-3680

Optional pump
Model GRZ-4610

Model:	GRZ™
Supports:	up to 1,000 lbs
Size:	36"W x 80"L x 7"H
	36"W x 76"L x 7"H also available in
	39"W, 42"W, 48"W, 54"W 60"W
Bladders:	9 self adjusting air cells
Fire Code:	16 CFR 1633
Cover:	Vyvox I™ Top, Medical grade vinyl bottom
Warranty:	5 years Limited
PDAC:	E0277 (requires alternating pump)

RAPID-AIR™ and RAPID-AIR PLUS™

Alternating Air Floatation Mattress Systems with Gentle Low Air Loss provides pressure redistribution for the prevention and treatment of Stages I-IV pressure sore ulcers. The model 4300P Rapid-Air Plus™ features a PPF™ 2.5" pressure redistribution foam base that protects the patient from bottoming out during power outage. Pediatric model 4300PM. Model 4300SRS with Side Rail Safety System.



Model 4300

Model:	4300 / 4300P / 4300PM
Supports:	275 lbs.
Size:	35"W x 80"L x 5.5"H (standard)
	25"W x 45"L x 3.5"H (pediatric)
Bladders:	18, 5.5" air cells
Fire Code:	Cal. Tech Bulletin #117
Cover:	***Vyvox III™
Warranty:	1 Year Limited
PDAC:	E0277


MATTRESS OPTIONS

SIDE RAIL SAFETY SYSTEM - SRS


SRS Safety Rail System can be added to all Alternating Pressure Mattresses for additional patient safety. To reduce accidental falls.

OPTIONAL INFLATOR

Blower pump can be used through the CPR Valve to inflate the mattress in under four minutes.

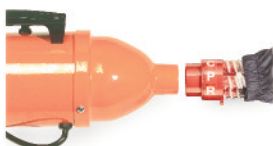

Model 9623
Optional to Inflate/Deflate

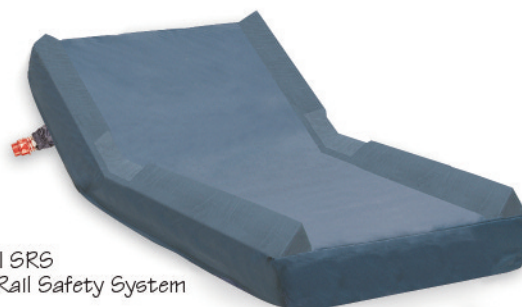

Model SRS
Side Rail Safety System


LATERAL ROTATION WITH LOW AIR LOSS

POWER TURN™ ELITE

The Power-Turn™ Elite Mattress System combines two clinically effective therapies: Lateral Rotation with Low Air Loss in a portable blower based system provides proper pressure redistribution. Lateral Rotation Therapy reduces the risk of pulmonary complications for patients who are immobile. A Full, Left or Right Turn mode is offered with either a 25° Small Turn or a 40° Big Turn.



Model 9800

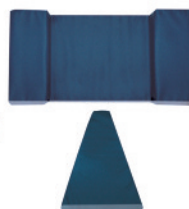

Optional Head Rest and
Leg Abductor Wedge

Model:	9800
Supports:	650 lbs. - 1,000 lbs.
Size:	36"W x 80"L x 9"H avail. in 42"W, 48"W
Bladders:	18 - 9" air cells, ventilated, low air loss
Fire Code:	Cal. Tech Bulletin #117
Cover:	***Vyvox III™
Warranty:	1 Year Limited
PDAC:	E0277

*Vyvox I™ is a soft non-shear nylon, fluid and stain resistant, anti-microbial and fire retardant.

THERAPEUTIC FOAM SLEEP SYSTEMS





SUSPENSION
Model BCSS II



SUPREME
Model BCSP II



ELITE
Model BCEL II

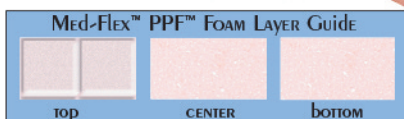
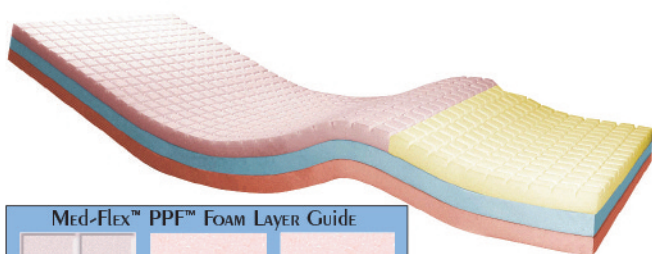


PRO
Model BCPR II

PRE-VENT II SERIES®

The Pre-Vent II Series® features five pressure redistribution foam mattresses. The Pre-Vent II Series® incorporates Blue Chip Medical's exclusive Sphere™, Reflect™ SCT™ and Multi-Flex™ foam contour designs. These react to body weight displacement across the entire surface to prevent and treat pressure sores. Blue Chip's unique construction together with our Med-Flex™ PPF™ foams, react instinctively to a compression set, allowing for a fully automatic, totally zoned pressure-redistribution surface. All models feature the exclusive Heel Wrap™ PPF™ foam heel section and can be ordered with a special upgraded fire barrier protection to meet flame retardancy specifications. ****This Includes California Technical Bulletin, #129, 603, Life Safety Code 101 and City of Boston 1X-11. Custom sizes and bariatric mattresses to support patients up to 1000 lbs available.

Model:	Suspension, Supreme, Elite, Pro, RC
Supports:	350 lbs, Bariatric 650 lbs - 1000 lbs. Elite only
Size:	36"W x 76"L x 6"H, (80"L & 84"L)
Addl. Widths	39"W x 42"W x 48"W (Elite Only)
Core:	Med-Flex™ PPF™ - Reflect™ SCT™
Fire Code:	16 CFR 1633****
Cover:	BCSS - Vyvex III™, BCSP - Vyvex II™ BCEL - Vyvex II™, BCPR - Vyvex I™, BCRC - Vyvex I™
Warranty:	BCSS 8 Yr., BCSP 7 Yr., non pro-rated BCEL 6 Yr., BCPR 6 Yr., BCRC 5 Yr., non pro-rated

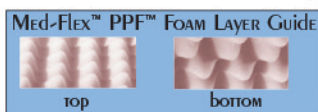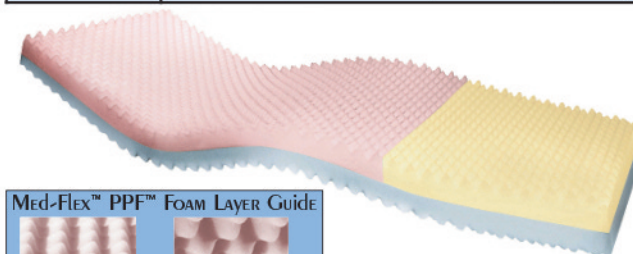

RC II
Model BCRC II

HEEL WRAP




Reflect



PRE-VENT II Shown In Cover



Perimeter Safety Border
(PSB) Included

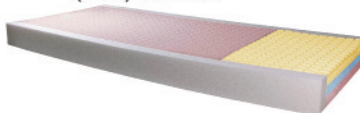

Shown with Optional
Safety Rail System (SRS)



RELIEF-CARE™ PRO

The Mattress System uses advanced Surface Cut Technology (SCT™) to create a clinically effective support surface. SCT™ enables Blue Chip to produce a therapeutic dual-zone PPF™ Pressurized Polymer Foam support surface with varying features and properties. Blue Chip's design absorbs and redistributes patient's weight to provide greater comfort, maximum support and proper pressure redistribution around bony prominences.



Model BCHTM

Model:	BCHTM
Supports:	250 lbs.
Size:	36"W x 76"L x 5.5"H (80"L, 84"L) 42"W x 80"L x 7"H
Core:	PPF™ Pressurized Polymer Foam
Cover:	Vyvex I™ Top, Medical grade vinyl bottom
Fire Code:	16 CFR 1633
Warranty:	1 Year Limited
PDAC:	E0184

THERAPEUTIC FOAM SLEEP SYSTEMS

HOME-CARE ELITE™

Designed for use in Patient Homes, Hospitals, and Institutions. The inner core is a high density, fire-retardant polyurethane foam that is covered with anti-microbial fire-retardant material. The mattress is closed with a self-ventilating zipper, and can be ordered in a variety of different sizes.



Model BCHM

Model:	BCHM
Supports:	250 lbs.
Size:	36"W x 76"L x 5.5"H (80"L, 84"L)
	42"W x 80"L x 7"H
Core:	High Density Polyurethane Foam
Cover:	Chemsafe™ Vinyl
Fire Code:	16 CFR 1633
Warranty:	1 Year Limited

PSYCHIATRIC AND PRISON MATTRESS

The inner core is a high density, fire-retardant polyurethane foam that is covered with heavy duty anti-microbial, fire-retardant vinyl stitched with super-strong nylon thread and has sealed seams. The mattress can be special ordered for use in prisons and psychiatric areas.



Model BC-PSY

Model:	BC-PSY
Supports:	250 lbs.
Size:	36"W x 76"L x 5.5"H (80"L, 84"L)
	42"W x 80"L x 7"H
Core:	High Density Polyurethane Foam
Cover:	Heavy Duty Chemsafe™ Vinyl
Fire Code:	16 CFR 1633
Warranty:	1 Year Limited

CRIB MATTRESS

Designed for use in Hospitals, Patient Homes and Institutions. The inner core is a high compression, fire-retardant polyurethane foam that is covered with Blue Chip's Chemsafe™ anti-microbial, fire-retardant cover. The mattress is closed with a self ventilating zipper, or under special order, a French Seam.



Model BC-PED

Model:	BC-PED
Size:	26"W x 50"L x 4"H, 27"W x 52"L
	27"W x 58"L, 32"W x 68"L
Core:	High Density Polyurethane Foam
Cover:	Chemsafe™ Vinyl
Fire Code:	16 CFR 1633
Warranty:	1 Year Limited

GEL-PRO™ SLEEP SYSTEM 5

The 3 section, 5.5" mattress utilizes a high density, high compression foam base designed to extend the life of the mattress and maximize its clinical benefit. 80% gel surface coverage. The zippered cover allows for sections to be removed for easy portability and storage. Custom sizes available. Bariatric sizes available by weight.



Model 6305

Model:	6305
Supports:	350 lbs
Size:	36"W x 80"L x 5.5"H
	42"W x 80"L, 48"W x 80"L
Core:	Gel / High Density Foam
Fire Code:	16 CFR 1633
Cover:	Vyvetex I™ Top, Medical grade vinyl bottom
Warranty:	1 Year Limited

MATTRESS OVERLAYS

AIR-PRO™ PLUS, AIR-PRO™ ELITE

The Air-Pro™ Plus & Elite Alternating Pump and Overlay Mattress Systems feature a heavy duty deluxe bubble style pad with head and foot tuck under flaps. Pads are available with or without end flaps and lasered air holes for ventilation. The adjustability of the Elite allows this system to be used for heavier patients. Multipad with high output 4400S pump supports patient weights to 350 lbs.



Model 4200
Model 4400

Model:	Plus 4200, Elite 4400
Pads:	4800, 4810, 4820, 4830
Supports:	250 lbs., 4830 supports 350 lbs.
Size:	35"W x 79"L x 2.5"H
Mattress:	Bubble style with flaps
Fire Code:	Cal. Tech Bulletin #117
Warranty:	Pump 1 Year, Mattress 6 months
PDAC:	E0181

GEL-PRO™ SLEEP MAT

The Gel-Pro™ mattress overlay offers maximum protection to guard against the incidence of pressure sores. Our viscous gel, migrates through sectioned bladders to provide proper immersion and envelopment of the bony prominences with over 80% surface coverage. A high density foam base supports the entire mattress overlay for maximum weight distribution and durability.



Model 6103

Model:	6103
Supports:	250 lbs.
Size:	35"W x 76"L x 3"H
Bariatric:	42", 48", 54" & 60" x 80" X 4" (based on weight)
Cover:	Vyvox I™ Top, Medical Grade vinyl bottom
Core:	Viscous Gel / Foam
Fire Code:	Cal. Tech Bulletin #117
Warranty:	1 Year Limited
PDAC:	E0185



Model 5500
Model 5700

STAT-H2O™, STAT-GEL™

The Stat-H2O™ water mattress helps to redistribute pressure while providing a cooler more comfortable sleeping surface. The internal baffle design balances the patient in the center of the mattress while eliminating the traditional water mattress "rocking" effect. Stat-H2O™ easily converts to STAT-GEL™, a full gel mattress overlay with the addition of our special polymer gel powder.

Model:	STAT-H2O™ 5500, STAT-GEL™ 5700
Supports:	250 lbs.
Size:	32"W x 72"L x 3.5"H
Bladders:	Medical Grade Vinyl with Baffles
Warranty:	6 Months
PDAC:	5500 E0198, 5700 E0185

STAT-AIR™

The mattress is constructed using a pliable medical grade vinyl that contours to the patient's body. Heat build up and moisture are greatly reduced, due to the air holes throughout the mattress. For patient safety adjustable anchor straps secure Stat-Air™ to an existing mattress. Quick release valve for rapid deflation during medical emergencies. The mattress is anti-microbial and fire retardant.



Model 5300

Model:	5300
Supports:	250 lbs.
Size:	35"W x 72"L x 3.5"H
Mattress:	Medical Grade Vinyl
Warranty:	6 Months
PDAC:	E0197

REPLACEMENT COVERS

Replacement covers extend the life of your mattress. They are available for all mattresses with removable Vyvox-I™, Vyvox-II™ or Vyvox-III™ top covers or Chemsafe™ medical grade vinyl. They are fluid proof and offer a low-shear, low-friction surface and effectively protect against a wide spectrum of micro-organisms, including MRSA, mold, mildew, staining and odors. These covers feature Blue Chip's exclusive Barrier Lok™ infection control system.



*Vyvox I™ is a soft non-shear nylon, fluid and stain resistant, anti-microbial and fire retardant.

**Vyvox II™ is 70 denier nylon, butyl backed, fluid-proof and stain resistant, anti-microbial, fire retardant.

***Vyvox III™ is multi-stretch knit urethane, low shear, fluid and stain proof, anti-microbial and fire retardant.